

Our Mission

Fostering a culture of innovation.
Leading healthcare solutions.
Simplifying complex care.



www.soleohealth.com

Administration of Elastomeric Device Medications




ADMINISTRATION OF ELASTOMERIC DEVICE MEDICATIONS

Your medication has been pre-mixed by your Soleo Health pharmacist and placed in an elastomeric device. Store all medication as directed. Some medications need to be refrigerated for stability.

- Check the label of each dose for your name, name of the medication, volume and the expiration of the drug.
- Each dose should be free of cracks, leaks, floating objects and discoloration. If you have any doubts call your Soleo Health pharmacist.
- Remove one dose from your refrigerator one (1) hour prior to infusing to allow the medication to warm to room temperature. **DO NOT HEAT OR MICROWAVE** the medication.
- Use each elastomeric device only once.

SCAN THE QR CODE BELOW TO VIEW THE VIDEO

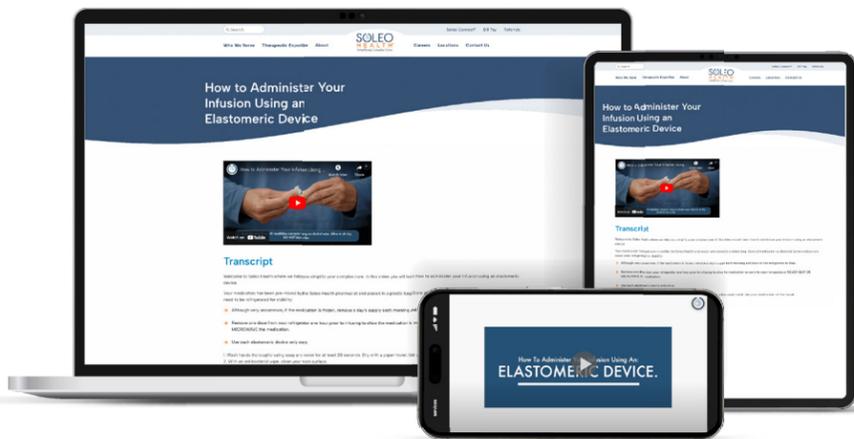

STEP 1

Wash hands thoroughly using soap and water for at least 20 seconds. Dry with a paper towel. Use your towel to turn off the faucet.

STEP 2

Clean your workplace and gather your supplies.

1. ELASTOMERIC DEVICE



2. ALCOHOL WIPES



3. FLUSH SYRINGES



STEP 3

Inspect your IV site for:

REDNESS

SWELLING

TENDERNESS

WARMTH

**IF THERE IS A CHANGE,
DO NOT INFUSE AND
CONTACT YOUR NURSE**

STEP 4

Inspect your medication.

IS IT EXPIRED?

IS YOUR NAME ON
THE LABEL?

IS IT THE RIGHT
MEDICATION?

ARE THERE ANY VISIBLE SPECKS OR FLOATING PARTICLES
IN THE MEDICATION?



STEP 5

Swab the end of your needleless connector using an alcohol wipe. Allow to air dry. Do not blot or fan.



STEP 10

Disconnect the tubing from your needleless connector.



STEP 6

Remove the cap from the saline flush. Holding the syringe upright, gently pull the plunger back to loosen it, and then push the plunger up to remove the air from the syringe. Attach saline flush syringe to catheter adapter by firmly pushing and twisting. Gently but firmly flush your catheter with ____ mL of saline.



STEP 11

Swab the end of your needleless connector using an alcohol wipe. Allow to air dry. Do not blot or fan.



STEP 7

Remove the protective cap from the end of the elastomeric tubing and attach to your needleless connector.



STEP 12

Flush the IV catheter with saline flush solution and swab the end of your needleless connector using an alcohol wipe as before. If instructed, follow with a Heparin flush using units/mL ____ mL. Call your nurse if you cannot flush your catheter or meet resistance.



STEP 8

Unclamp the pinch clamp on the elastomeric tubing to start the infusion.



STEP 13

Discard the used elastomeric device in the appropriate waste container.



STEP 9

When the infusion is complete, close the pinch clamp on the elastomeric tubing before disconnecting.



